


DRUG STABILITY STUDY:TOBRAMYCIN

Low Level

TOBRAMYCIN-mg/L			
Time (h)	Vacurette Serum	Vacurette Serum Gel	Vacutainer Plain Red
0	0.5	0.5	0.5
4	0.5	0.5	0.5
24	0.4	0.5	0.4
48	0.6	0.5	0.5
Mean	0.5	0.5	0.5
CV (%)	11.7	8.0	7.1
SD	0.1	0.0	0.0

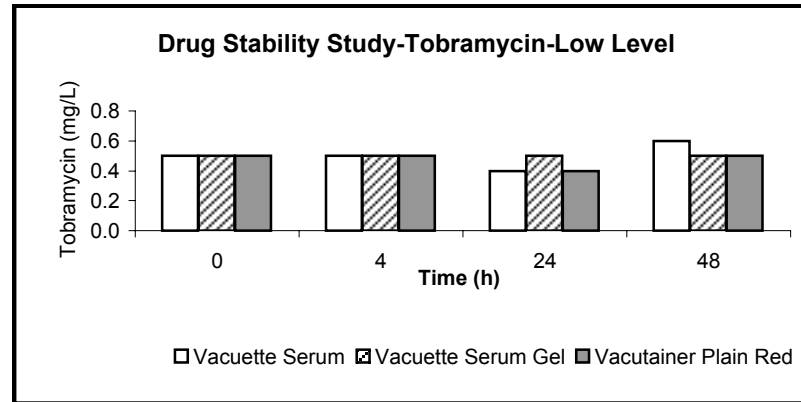

High Level

TOBRAMYCIN-mg/L			
Time (h)	Vacurette Serum	Vacurette Serum Gel	Vacutainer Plain Red
0	5.1	5.0	4.8
4	5.2	5.2	5.1
24	4.7	4.4	4.3
48	4.6	4.4	4.1
Mean	4.9	4.8	4.6
CV (%)	6.2	8.1	10.2
SD	0.3	0.4	0.5

