


DRUG STABILITY STUDY: VALPROIC ACID

Low Level

VALPROIC ACID-umol/L			
Time (h)	Vacurette Serum	Vacurette Serum Gel	Vacutainer Plain Red
0	152	157	154
4	152	154	149
24	152	159	152
48	152	157	160
Mean	152	157	153
CV (%)	0	1	3
SD	0	2	5

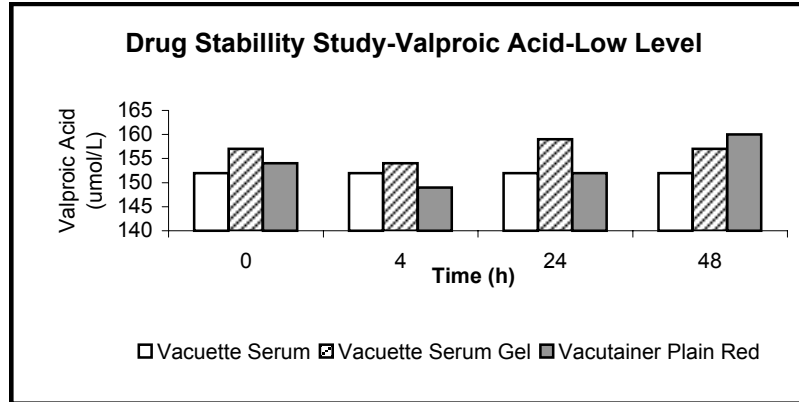

High Level

VALPROIC ACID-umol/L			
Time (h)	Vacurette Serum	Vacurette Serum Gel	Vacutainer Plain Red
0	504	496	475
4	534	505	499
24	500	513	504
48	509	522	501
Mean	512	509	495
CV (%)	3	2	3
SD	16	11	13

