


## DRUG STABILITY STUDY:VANCOMYCIN

### Low Level

VANCOMYCIN-mg/L			
Time (h)	Vacurette Serum	Vacurette Serum Gel	Vacutainer Plain Red
0	8.2	8.0	8.5
4	8.4	7.8	8.5
24	7.5	8.0	8.1
48	7.9	7.6	7.5
Mean	8.0	7.9	8.2
CV (%)	5.3	2.2	5.9
SD	0.4	0.2	0.5

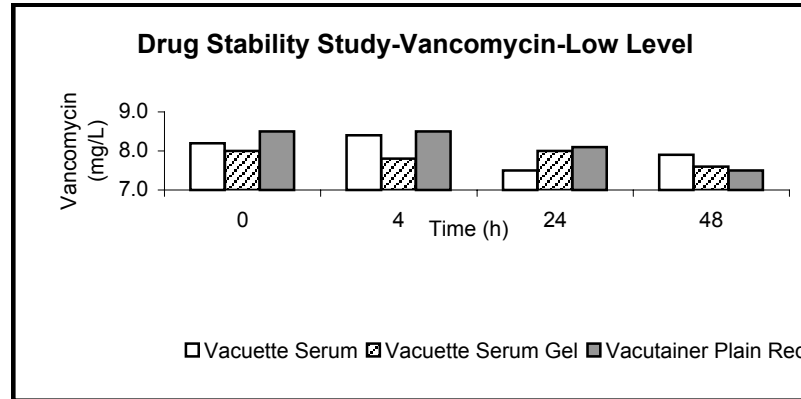

### High Level

VANCOMYCIN-mg/L			
Time (h)	Vacurette Serum	Vacurette Serum Gel	Vacutainer Plain Red
0	33.7	34.5	33.4
4	35.0	33.7	34.2
24	30.2	33.7	34.8
48	33.7	33.4	33.5
Mean	33.1	33.8	34.0
CV (%)	6.2	1.4	1.9
SD	2.1	0.5	0.7

