

3 PART SYRINGES- Luer Lock/ Slip/ Catheter Tip

Reference Items	Description
REF:	SY3-XLC/SC/CTC-EC
Associated group	Injection Systems
Product description	Syringe
Intended use	Injection therapy

Reference code	Description	Size (ml)	Type	Scaling (ml)
SY3-50ESC-EC	3 part syringe 50 ml Luer Slip	50-60ml	ECCENTRIC	1
SY3-50LC-EC	3 part syringe 50 ml Luer Lock		CENTRIC	
SY3-50CTC-EC	3 part syringe 50 ml Catheter Tip		CENTRIC	




Luer Slip – 50ml-
Graduation until 60ml



Luer Lock – 50ml-
Graduation until 60ml



Catheter tip- 50ml –
Graduation until 60ml

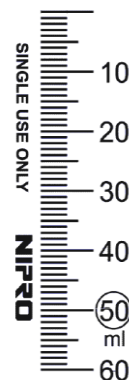



Barrel - Scale

PRODUCT DATA SHEET

<u>Product Registrations</u>	
CE- marked	CE
Classification	Medical Devices
Medical Product class	Class I
Rule	Rule 2
Guideline	MDD 93/ 42/ EEC
UMDNS	13929
<u>Standard norms:</u>	
ISO 7886-1/ ISO 15223-1/ ISO 11607-1/ ISO 11607-2/ EN ISO 9001 EN ISO 13485/ ISO 10993-1/ ISO 10993-4/ ISO 10993-4/ ISO 10993-5/ ISO 10993-7/ EN 1707/ ISO 594-1/ ISO 594-2	


<u>Manufacturing details</u>	
Legal Manufacturer	Nipro Europe
Country of origin	China

<u>Product materials :</u>	
All materials used to produce our 3 part syringes are medical grade, latex-free and PVC-free.	
Polypropylene	-Plunger, Barrel - Hub and protector of the Catheter tip
Isoprene rubber	Gasket
Silicone	Lubricant
Epoxy	Glue


<u>Product Sterilization</u>		
Sterile	ETO Sterilization procedure Ethylene – Oxide	
Do NOT re-use	Single use only	
<i>Each outer box (or shipping carton) contains a chemical indicator, which ensures sterility. This indicator is a blue sticker, which turns red when sterilized.</i>		

PRODUCT DATA SHEET

Packaging

- LOT-number
- Expiry : 5 years 
- REF: SY3-50SC/LC/CTC-EC
- QTY- pieces
- Barcode

Example of Label:

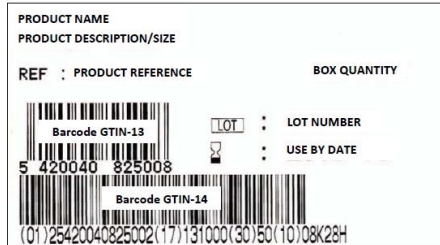

An individual package, is sterile paper for blister package. Inner box is made of duplex board, and outer box of double corrugated board.

Transport conditions

Closed and dry

Storage conditions

Only open the packaging before use to guarantee sterility.

Reference	Size (ml)	Scaling	Quantity inner box	Quantity Outer box	Net Weight (g)	Gross Weight (kg)	INNER BOX			OUTER BOX			
							L mm	W mm	H mm	L mm	W mm	H mm	M ³
SY3-50SC-EC	50	1	30	360	47,000	17,000	280	210	140	580	440	460	0,11
SY3-50LC-EC	50	1	30	360	47,000	17,000	280	210	140	580	440	460	0,11
SY3-50CTC-EC	50	1	30	360	47,000	17,000	280	210	140	580	440	460	0,11

15 Languages on packaging according the specifications of European Medical Device Directive:

- English, French, Dutch, German, Spanish, Italian, Portuguese, Greek, Swedish, Danish, Norwegian, Finnish, Russian, Turkish, Polish

**Contact details**

Nipro Europe NV
Weihoek 3H
1930 Zaventem
Belgium

Phone: +32 (0)2 725 55 33 / Fax: +32 (0)2 725 70 41

E-mail: info@nipro-mediceurope.comWeb: www.nipro-mediceurope.com